



## Baby Incubator BB-200(Standard)

### Function features

- Microprocessor based servo controlled temperature system
- Control modes: Air mode and baby skin mode
- Humidity is adjustable continuously
- Set temperature, air temperature, baby skin temperature, timer and heating power are displayed separately by LED
- Integral timer to record therapy time
- Self - testing function, various failure alarms by audible and visual
- Alarm function:  
Power failure, Over temperature, Temperature deviation, Fan failure, Temperature sensor failure
- >37° temperature set function
- Triple protection for over temperature with separate cut off device, more safety system
- The inclination of infant bed is adjustable
- Single wall hood, 4 operating windows and 2 iris ports, 5 infusion seals
- Independent locking device for front door
- RS232 connector, Oxygen inlet, 1 tray

## Technical data

Power supply	AC110/220V 60/50Hz
Power input	650VA
Environment temperature	20°C ~ 30°C
Operating condition Environment relative humidity	30% ~ 75%
Environment air velocity of flow	< 0.3m/s
Air temperature control range	25°C ~ 37°C 37.1°C ~ 38°C
Skin temperature control range	32°C ~ 37°C 37.1°C ~ 38°C
Skin temperature sensor accuracy	±0.3°C
Temperature fluctuation	±0.5°C
Uniformity of mattress temperature	≤0.8°C
Warm-up time(from25°C)	< 30min
Internal noise level	< 50dB(A)
Infant bed tilt angle	±10°
Mattress size	65cm(L)*37cm(W)
Trough capacity	1200mL
Air filter	0.5μm

## Transport and storage

### Package

Each unit is packed in two plywood cases  
(control unit and cabinet):112\*66\*98cm

N.W.:80Kg G.W.:110Kg

Hood:98\*63\*74cm

N.W.:15Kg G.W.:30Kg

### Transport and storage

Environment temperature -10°C ~ +70°C

Environment relative humidity ≤80%

Atmospheric Pressure 500 ~ 1060hpa

## Standard configuration

Main body(including the Transparent hood,Control system,Infant bed,Water trough),I.V.pole,Mattress, Drawers,Castors,Tray,Oxygen inlet



## Optional functions




LCD Touch Screen



X ray cassette tray&Integrated scale



Electronic Elevator

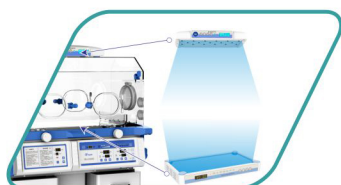

Upside LED phototherapy unit BL-30D;  
Downside LED phototherapy unit BL-100D;



Double wall hood



Air&Oxygen Blender

ZHENGZHOU DISON INSTRUMENT AND METER CO.,LTD.

Add: No.67 Jinsuo Road, High-Tech Development zone, Zhengzhou, 450001, Henan, China

Tel: 0086-371-67897160

Fax: 0086-371-67897164

E-mail:info@zzdison.com

**disonmed**  
www.zzdison.com